User Manual

FY-6


2-in-1 Shiatsu Foot and Back Massager with Heat

Thank you for purchasing 2-in-1 Shiatsu Foot and Back Massager with Heat.
With normal care and proper treatment, it will provide years of reliable service.
Please read all instructions carefully before using this product.
Please retain this manual for future reference.

FEATURES

- Shiatsu, rolling massage on sole of the foot
- With optional heat function
- 3 massage levels setting
- 3 massage modes
- Detachable and washable cover


CONTENTS

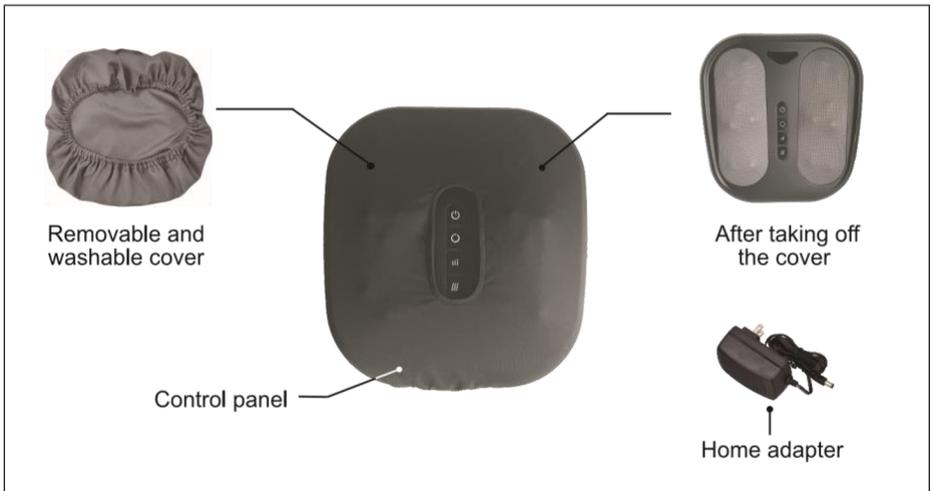TECHNICAL DATA

Measurements:	13.4 x 14.4 x 3.55 inches
Weight:	4.4lbs
Voltage:	Input: AC 100-240V ~50/60Hz Output: 12VDC 2A
Nominal Power:	max. 24Watt
Automatic runtime:	15 minutes

1. Please put the product on the flat floor.
2. Connect the adapter cable to the corresponding port on the foot massager.

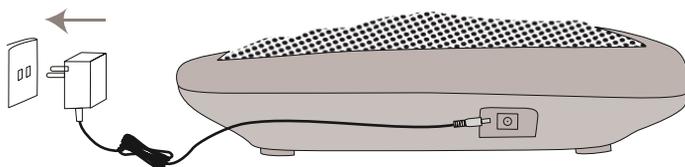3. Plug the home adapter into electrical outlet (the adapter is safe to use with both 110-120V & 220-240V)
4. Turn on the massage device through operating control panel (instructions on page 7).

- This device will switch off automatically after 15 minutes timer run out.
- Do not setup or use the device in a bathroom or similiary wet / damp areas.

- Before plug ,please check and make sure the power cord is undamaged.
- Do not use the wet hand to plug in the adapter.
- Please don't use the product by standing posture or put any heavy object on the product to avoid any broken caused.
- Please don't stand or seat on the product .
- When moving the product, please don't pull the power cord or drag the product directly to avoid any damaged caused.
- Please don't let children and pets play the products for preventing any accident happened.
- Do not use the wet hand to pull out the power cord in order to avoid any electric shock.
- Do not use this product in damp environment to avoid electric shock.

Please read the following instructions carefully before putting your massage device into use to ensure trouble-free operation and optimum efficiency. Please retain these instructions for further reference!

- The device conforms to recognized technical principles and the latest safety regulations.
- This item is NOT A TOY. Close supervision is necessary when this appliance is used by, on, or near children or disabled persons. Do not allow children to play with product.
- Do not allow more than one person to use this product at a time.
- The device should never be left unattended when plugged in.
- Never use any power sources other than UL standard home adapter provided as original equipment with this massage device.
- Any possible repairs may only be carried out by authorized expert staff. Improper use and unauthorized repairs are not permitted for safety reasons and lead to loss of warranty.
- Please avoid contact of the device with water, high temperatures and direct sunlight.
- Do not damage, process, excessively bend, pull hard, twist or knot power cord.
- Do not use any damaged cables, plugs or loose sockets.
- Do not use the product when the power cord is damaged.
- In case of malfunction, disconnect immediately from the mains.
- Misuse or incorrect use excludes any liability for damage.
- Do not use while you sleep.
- In order to avoid excessively stimulating the muscles and nerves, the recommended massage time should not be exceeded 15 minutes at a time.
- Do not use this product for medical treatment.
- People who are insensitive to heat must be accompanied by normal person when using this massage device.
- During use, if you feel abnormal, please stop using this product immediately, and then consult your doctor.
- Do not use this product for medical treatment.

We recommend that those with any of the following conditions consult a doctor before using this massage device.

1. Those who are using electronic medical instrument imbedded in the body, such as cardiac, pacemaker etc.
2. Those who are being treated by doctors, especially those who feel uncomfortable.
3. Those who are malignant tumor patients, heart disease patients, acute patients.
4. Women who are pregnant or menstruating.
5. Those with osteoporosis or fracture spine.
6. Those with skin disease or those whose skins are injured.
7. Those whose body temperature is over 38 degrees (Febrile Phase).

Failure to observe the above instructions may constitute misuse of the product and can cause serious injury or burns.

- Please choose a chair with desired height and be seated, then put your feet on the product (see right picture). Operate below control panel to get foot massages.

**Power Button**

Press to turn on the massager, the indicator turns Yellow.
The massager activate the auto mode; Press to turn off the massager.

**Direction Button**

Direction Button When the massager is turned on, Press to shift Shiatsu clockwise rotation(Green indicator) or Shiatsu anti-clockwise rotation(Red indicator)

**Speed button**

Low Speed --Indicator turns Green.
Medium Speed--Indicator turns Yellow.
High Speed--Indicator turns Red.

**Heat Button**

High Level - Indicator turns to highlight red
Low Level - Indicator turns to lowlight red.

After use

- Turn off the power button on the control panel.
- Remove the plug from the power supply socket.

Caution

- Please don't put water or other liquid on the control panel.
- Please don't press the button too long time at a time to avoid the program confused.
- When power failure happens during massage, please turn off the power and pull out of the plug from socket.
- Please do not let your hand or fingers close to and inserted into the gap of products for preventing any injuries caused.
- Please do not let children and pets play the products for preventing any accident happened.

- Before cleaning, make sure the power supply is turned off and the plug has been pulled out from the socket.
- It is forbidden to dip the whole product in water or other liquid when cleaning.
- Use a soft cloth to wipe the surface of the product and control panel.
- Use a piece of cloth that has been dipped with cleaning agents to wipe the main unit of the foot Massager, and then use a piece of dry and soft cloth to wipe it.
- Do not pull the power plug by damp hands.
- Do not use water washing product, do not spray the water or other liquids to the interior of the product.
- Do not use corrosive cleaning agents like benzene or thinner to clean the Foot Massager, otherwise, the color of upholstery will fade or the Foot Massager will crack or be eroded.
- Do not iron the upholstery.
- When product is wet, let the product dry naturally and thoroughly.

Removable & Washable Cover

- Detach the cover and wash in machine or by hand. Allow to air dry

Maintenance

- Make sure the power supply is turned off and the plug has been pulled out from the socket.
- Sweep off the dust. Keep the foot massager clean and away from damp places.
- Store the massager in a cool, dry place out of reach of children.
- Avoid contact with sharp edges or pointed objects which might cut or puncture the cover of massager.

TROUBLESHOOTING GUIDE

Fault	Solutions
There is a clattering noise in the movement	This is perfectly normal caused by the working of product. It is normal status.
The massage device switches off suddenly	Device may have been in use for 15minutes. Auto-shut off function activated. Continue using by pressing the 'power' button after resting the unit for 15minutes.
The massage device fails to start	Check if device has been used for more than 15mins. If so, let device rest for 15mins before continuing. Device may require longer time to cool in warmer environments.
	Check if adapter is plugged into socket.
	Check if adapter output is connected to device power inlet.

12 months Warranty

If the product break down within 12 months in the period of proper use, customers can enjoy the relevant warranty to get them replaced.

Attention!

No warranty will be given to any force majeure and man-made causes to a defective product, such as improper care, personal tear down and willful damage, etc.

FCC STATEMENT

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio / TV technician for help.

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

WARNING: Changes or modifications to this unit is not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

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