


LED FACIAL MASK

Important notes:

1. Please use the product according to the correct method; it is always better to read the user manual before using the device, so as to avoid malfunctions caused by incorrect operation.
2. Before wearing the LED mask, we advise you to carry out a light sensitivity test: simply switch the device on and place it on your arm or leg for 3 to 5 minutes; if your skin turns red or you feel any discomfort, stop using it and contact us for after-sales service.
3. Do not use the device if you are pregnant, breastfeeding or taking medication. Please consult your doctor before use.
4. Do not let children use it — it is not a toy.
5. Do not use it if you suffer from physical, sensory or mental disorders, or if you are undergoing dermatological treatment, including during a flare-up (such as lupus erythematosus, photosensitive eczema or albinism, etc.)
6. This product is not waterproof; do not immerse it in water, you can clean it with a damp wipe.
7. It is recommended to wear the protective goggles if you are sensitive to light. If you still feel any discomfort, we suggest closing your eyes and enjoying the session for 10 to 15 minutes.
8. It is recommended to use it 3 to 4 times a week, for 10 to 30 minutes per session.

Product operation explained:

- 1. Power button: press and hold for 1 to 2 s to switch on / off;
- 2. The E key is the energy adjustment key: press briefly to switch between the 3 power levels;
- 3. Timer key: press the T+ key to increase the lighting time
Press the T- key to reduce the lighting time
5 levels (10 / 15 / 20 / 25 / 30 min)
a single session of 15 to 25 minutes is recommended.
- 4. While charging, the battery indicator flashes; it goes out once charging is complete.

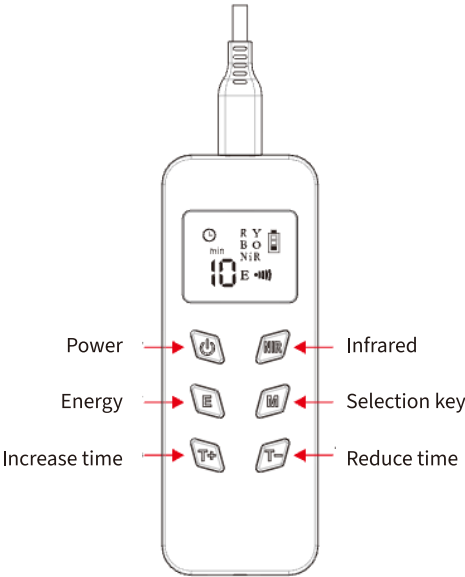

Four colours

			
Red light 630 nm	Blue light 460 nm	Yellow light 590 nm	Near-infrared light 850 nm

Packing list



Storage bag



Strap x3



Charging cable



Manual



The controller



Facial mask




Eye mask

Product specification:

Product name	LED facial mask device		
LED lamp beads	240	Carton size	34x23x4,5 cm
Timing	10/15/20/25/30 min	Frequency	0,3 A 50/60 HZ
Service temperature	0-40 C	Product net weight	0,32 kg
Adapter input	100-240 V	Adapter output	5V DC 3000MA
Four wavelengths of light	Red light 630 nm, blue light 460 nm, yellow light 590 nm and near-infrared light 850 nm		
Individual product list	Mask, controller, charging cable, manual, fixing strap		

How to put on the mask


		
Step 1	Step 2	Step 3

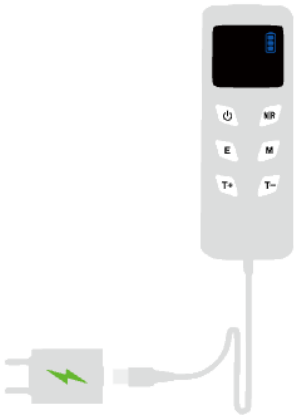

Step 4

Take out the controller and the USB cable. Connect them to each other. Plug the USB end into an adapter, then into a power socket.

Switch on the socket and charge the controller until 3 full bars are shown, as illustrated below:

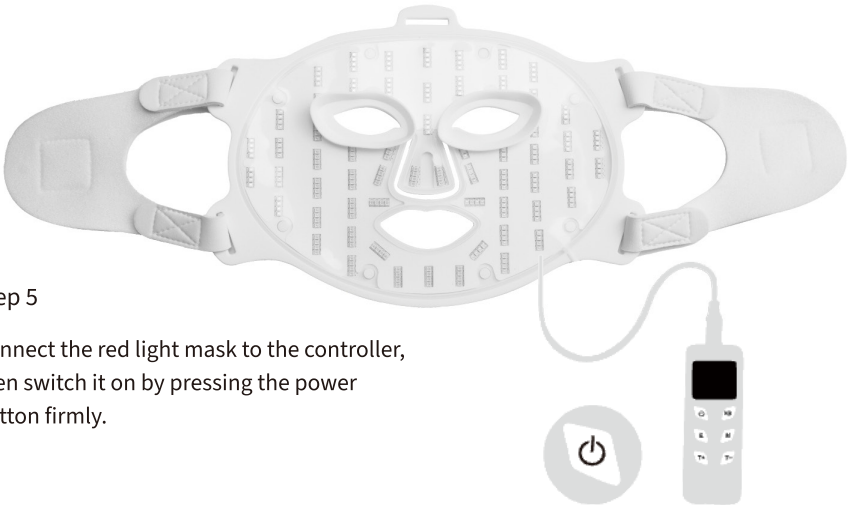

Step 5

Connect the red light mask to the controller, then switch it on by pressing the power button firmly.

LED FACIAL MASK:

The LED FACIAL MASK is a rechargeable LED light therapy device, worn at home for beauty care. The surface of this product is made of medical-grade silicone, convenient and easy to use, foldable and portable; it fits different faces thanks to the adjustable nylon buckle, which makes it suitable for different people. It consists of 240 high-quality LED diodes, divided into 60 blocks each containing 4 LED beads, with the four beads emitting 4 wavelengths (red light 630 nm, blue light 460 nm, yellow light 590 nm, near-infrared light 850 nm). The controller also allows 3 energy levels and different lighting times to be set. Based on light therapy, this product penetrates deeper into the skin, stimulates the collagen- and elastin-producing cells and achieves a beautifying effect on the skin.

Mode effectiveness:

Red light: 630 nm, anti-ageing, boosts collagen, firms, increases elasticity

Blue light: 460 nm, kills bacteria and acne bacteria, reduces inflammation and scarring, increases skin oxygenation

Yellow light: 590 nm, evens out skin texture, stimulates the red blood cells to reduce redness and fine lines

Near-infrared light: 850 nm, repairs the skin, improves blood circulation, increases skin hydration

The modes are: Red light, 1.2.3 Level

Blue light, 1.2.3 Level

Yellow light, 1.2.3 Level

The infrared light is controlled independently and can also be used with any other light.

First setting:

Before opening the packaging, always fully charge the mask battery; the battery indicator on the controller starts to flash to show that it is charging.

The battery indicator goes out once charging is complete.

Charge for 4 to 5 hours, then wear the mask comfortably using the adjustable velcro straps. To fit the mask, pass the strap through the opening on the other side of the mask, around the head, and secure it.