INSTRUCTION MANUAL

Red light neck device

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

Important: Change or modifications not expressly approved by the manufacturer responsible for compliance could void the user's authority to operate the equipment.

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
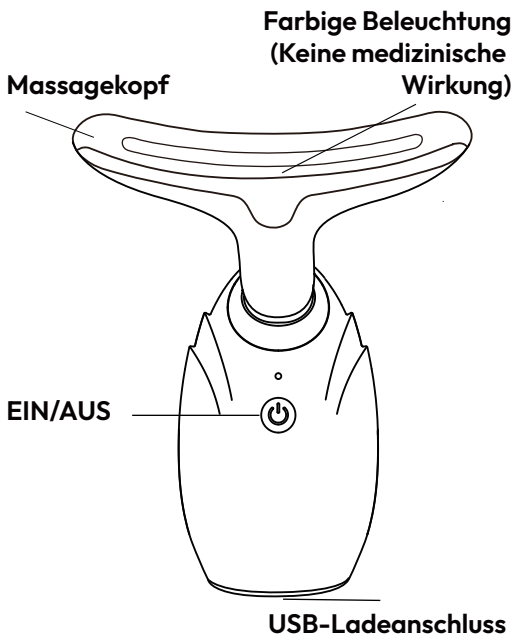
INTRODUCTION

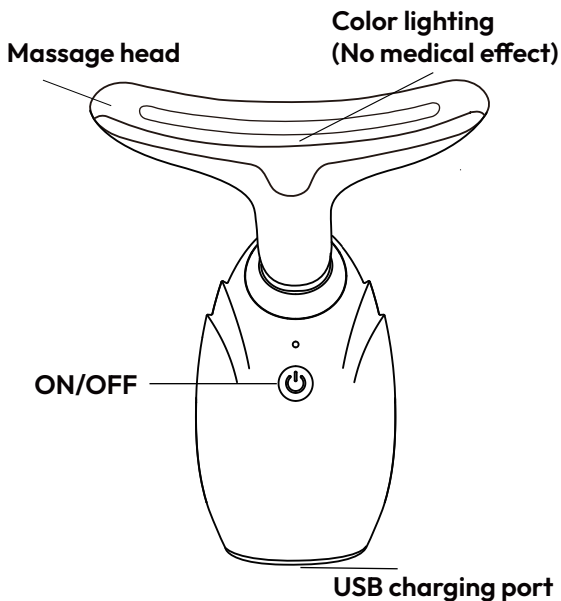
USAGE

1. After washing your face, apply your favorite skincare products evenly to your face.
2. Press the ON/OFF button for 2 seconds to turn on the device.
3. Short-press the ON/OFF button again to adjust your preferred lighting mode and massage mode
4. After using the product, press the ON/OFF button for 2 seconds to turn off the device.

WARNING

1. Not suitable for use during facial plastic surgery, subcutaneous metal implants, pregnancy, and menstruation.
2. Please wipe the skin rejuvenator with a tissue. Do not wash with water or put it together with metal objects.
3. Do not place it in a humid place such as the bathroom to avoid short circuits.
4. It can only be used on the face and neck, and it is not recommended for people with sensitive skin to use it every day but once a week.
5. When maintenance is required, please confirm whether the power is turned off first. If the power is turned on and a wet tissue or towel is used for maintenance, it may cause irritation.
6. Do not wash the body of the product; please wipe it with a wet tissue or a soft cloth dipped in a small amount of soapy water. Do not wipe the product with corrosive solutions (alcohol, banana oil, gasoline, acrylic remover, etc.).

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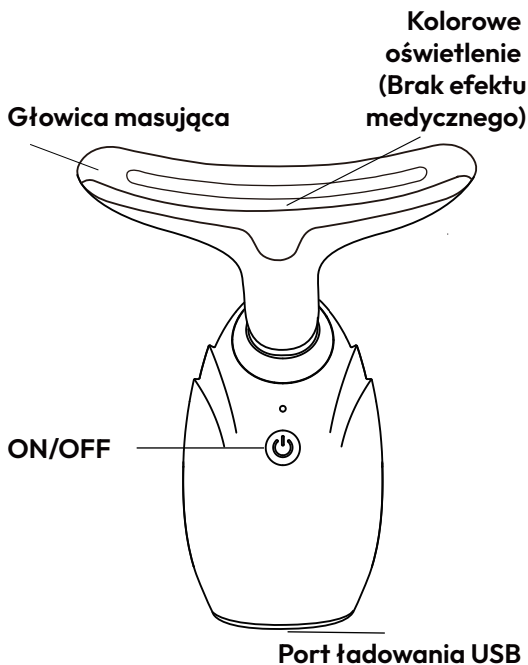DECLARATION OF CONFORMITY

Confinity NV, Dorp 16, 9830, Sint-Martens-Latem, Belgium
declares the following device in sole responsibility:

Brand name: Wellhealth
Product type: Red light neck device
Item number: OP_017013

Complies with the following harmonization rules:

EN IEC 55014-1:2021
EN IEC 55014-2:2021
EN 61000-3-3:2013+A2:2021
EN IEC 61000-3-2:2019+A1:2021

Signed for and on behalf of:

Sint-Martens-Latem, Belgium - March 2025
A. Pappijn - Product Manager