
## INSTRUCTION MANUAL

# Shiatsu Neck Massager

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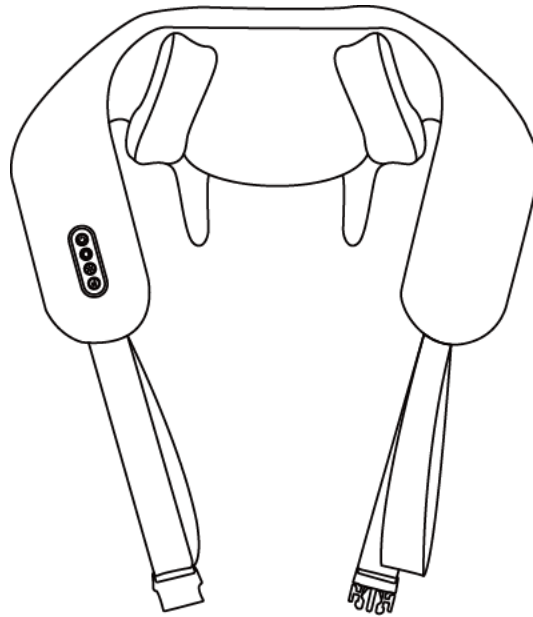## PRECAUTIONS

1. Please read the instructions carefully and do not use this product continuously for a long time to prevent muscle pain caused by excessive massage.
2. Do not use it when lying flat or half lying down.
3. Do not use this product during precise operation or driving.
4. Do not use if you have undergone surgery on the massaged part or have skin damage, redness, inflammation, congestion, and other symptoms.
5. People with heart disease or medical history, acute disease patients, skin disease patients, malignant tumor patients, bleeding tendency and other blood diseases, osteoporosis, soft tissue diseases, stroke, and other physical abnormalities should not be used.
6. Do not use it if you have slow movement, physical and mental discomfort, or intellectual disability. If you need to use it, please follow the doctor's advice.
7. Pregnant women and menstrual period are not recommended to use the massager. If you need to use it, please follow the doctor's advice.
8. When storing the product, please pay attention to placing it in a place where children are not easy to reach.
9. Children under the age of 8 are prohibited from using it, and children aged 8-12 must use it under the guidance and care of their parents or someone with supervisory ability.
10. . Do not scratch or stab this product with sharp objects.
11. This appliance is not intended for use by persons (including children) with reduced physical, sensory, or mental capabilities, or lack of experience and knowledge unless they have been given supervision or instruction concerning the use of the appliance by a person responsible for their safety.
12. Children should be supervised to ensure that they do not play with the appliance.
13. . It must only be supplied at safety extra low voltage corresponding to the marking on the appliance.


**Caution**

1. This product is not a medical device and is not suitable for diagnosis, treatment, and other commercial purposes.
2. Before use, please remove the jewelry on the massaged part or other items that may affect the massage.
3. During use, if there is pain or discomfort in the massage area, please stop using it.
4. When using, if the product is abnormal, please shut down immediately, cut off the power supply and contact the customer service hotline, do not disassemble or repair by yourself.
5. The massage area has a heating surface, people who are not sensitive to heat and need to pay attention when using it.
6. When using in other specific environments, please comply with safety regulations.

**Product Protection**

1. Please use the charging cable configured by our company to charge the product with a standard adapter, and pay attention to the safety of electricity use.
2. Please keep this product away from water, fire and corrosive environment.

**PRODUCT OPERATION**

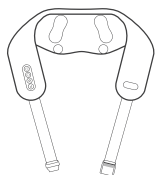Massager

User Manual

Charging Cable

**Necessary Knowledge**

**About frequency of use**

1. Use as needed, within 10 minutes at a time. It is recommended to take a proper rest before use after one use.
2. If overused, it may cause physical discomfort and skin injury.
3. Therefore, when you use it for the first time, you need to observe your physical condition and skin changes and shorten the use time as much as possible.

**Cautions Before Use**

1. When using this product, please take a personally comfortable posture and completely relax your body.
2. If you feel pain or discomfort on the skin, please stop using it and seek medical attention immediately.

**State After Use**

After long-term massage, the massage area may feel red and hot while feeling comfortable, which is normal.

**Cleaning**

If you need to clean the appearance, you can use a rag dipped in water or 75% alcohol to wipe the smooth surface, do not use water to rinse.

**FEATURES**

Bionic Massage




Free Hands

Gentle Heating

PRODUCT OPERATION

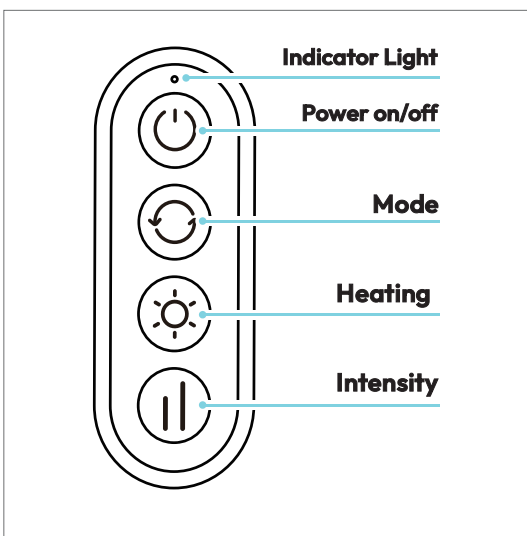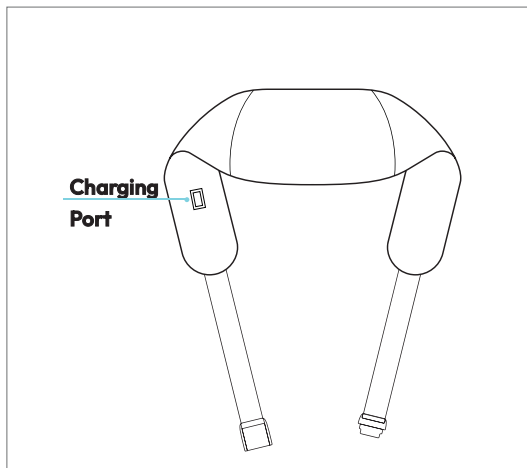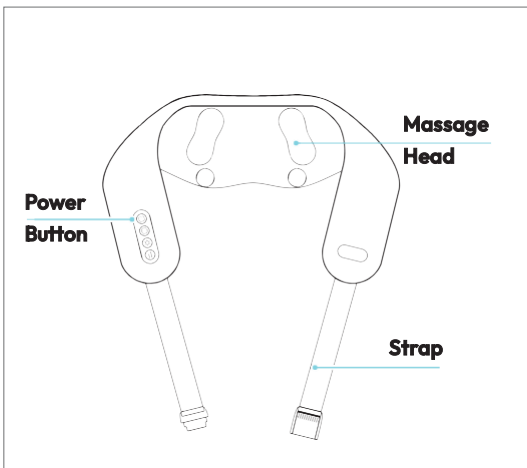**Correct Way to use**

**prepare before use**

Please check whether the product and accessories are complete and free of damage  
 Long press the power button for 2 seconds, and the product starts to work.

• **Mode Selection**

The default mode 1 is the default mode when the device is turned on. Short press the mode button to switch between mode 1 and mode 2

• **Heating**

The default heating is high-grade when the device is turned on, Short press the heating button to switch gears.  
 The sequence is high heating low heating turn off


• **Intensity**

The default intensity is the low intensity when starting up, short press the intensity button to switch between the low intensity and the high-intensity



• **Power Off**

Automatic shutdown: After 10 minutes of cumulative use, the product Stops working.  
 Manual shutdown: press and hold the power button for 2 seconds, the product will stop working.

## Charging

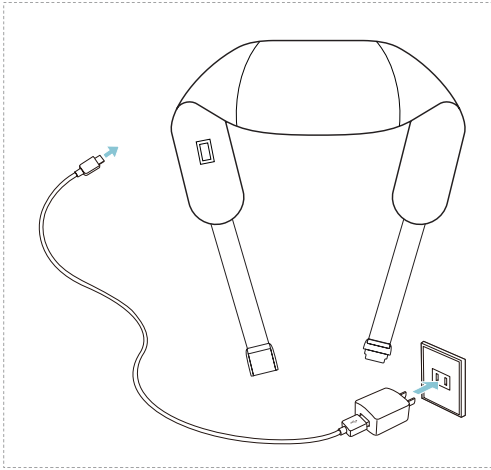1. This product cannot be used while charging
2. Insert the plug of the charging cable into the charging port on the port on the
3. When the product is working, it is connected to the power supply for charging, and the product starts a self-protection program and stops working automatically.
4. When the battery power is lower than 20%, the indicator light will show a steady red light.
5. After charging is complete, the indicator light will flash white. It is normal for the device and the power adapter to heat up slightly during charging.
6. After charging, please unplug the power.
7. This appliance contains non-replaceable batteries.
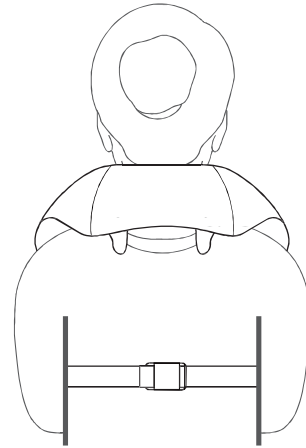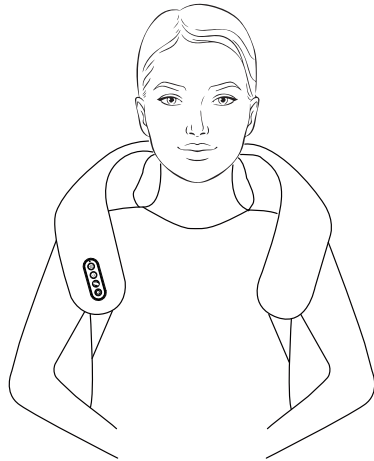
## Product Maintenance

If there is a failure during operation, please check the following. If the fault cannot be rectified or the product still cannot operate satisfactorily, please cut off the power supply and have it repaired by professional technicians. Do not disassemble the product by yourself!

## Disposal and Recycling

This product should not be thrown away with household waste, When you decide to dispose of this product, please follow local environmental laws and guidelines.

### Notice

- Before disposing of this product, the battery must be removed as required by regulations;
- When removing the battery, the product must be powered off;
- Batteries should be disposed of safely in accordance with laws and regulations.
- Non-professionals do not disassemble by themselves!

### Steps to Remove the Battery

- Cut off the external power supply and remove the product shell;
- Unplug the power supply on the battery;
- Disassemble the battery;
- Please place waste batteries at a fixed recycling point;
- Please abide by the local environmental protection laws and disposal guidelines, and care for the environment.

## Product Specifications

Product Name	Shiatsu neck massager
Model	OP_013095
Rated Voltage	5V <sup>==</sup>
Rated Power	25 W
Power Supply	7.4V <sup>=</sup> 2000mAh rechargeable lithium battery
Charging Time	≤ 3.5 hours
Working Time	≥ 60 mins
Executive Standards	GB4706.1-2005V GB4706.10-2008

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