

INSTRUCTION MANUAL

Water flosser

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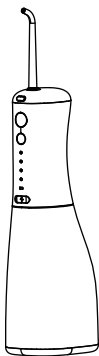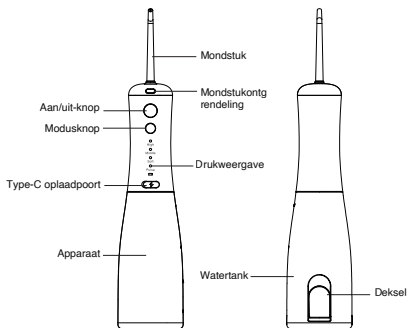Specificaties


Productnaam: Water Flosser

Model: F500

Spanning: DC5V, 1A

Batterij: 1500mAh lithiumbatterij

Waterreservoir capaciteit: 350 ml

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
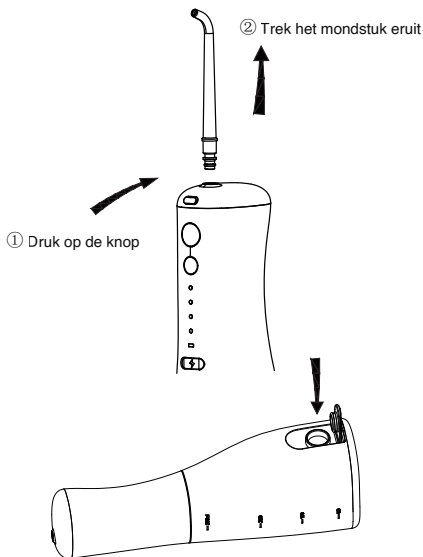
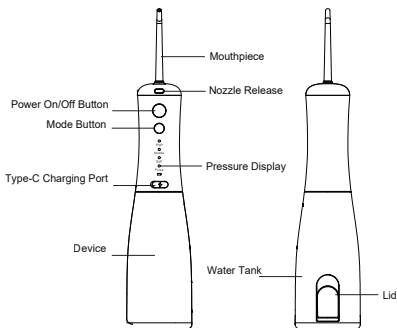
Battery: 1500mAh Lithium Battery

Water Tank Capacity: 350 ml

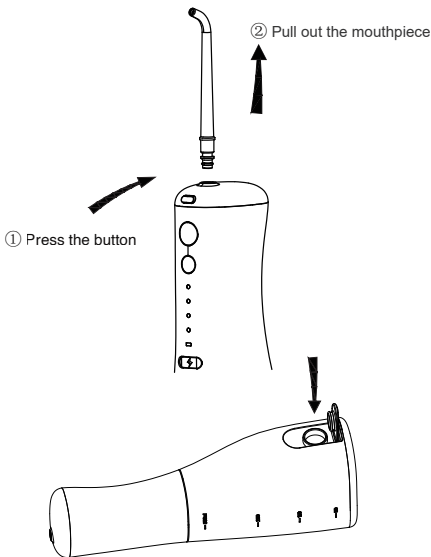1. Insert the nozzle into the hole at the top of the dental irrigator.
To replace the nozzle, first turn off the device and press the button for easy replacement.
2. Add water to the water tank:
Open the tank lid, hold the toothbrush horizontally, pour water into the tank, and then close the lid.
3. Lightly press the mode button to select your desired mode.
This product offers four modes with specific functions:

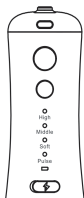- **High Mode:** The first mode indicator lights up, indicating strong water pressure for efficient cleaning. Suitable for most users.
- **Medium Mode:** The second mode indicator lights up, indicating moderate water pressure.
- **Soft Mode:** The third mode indicator lights up. Ideal for first-time users or those with oral sensitivity.
- **Pulse Mode:** The fourth mode indicator lights up, providing deep cleaning and promoting gum circulation.


4. Before use, hold the toothbrush vertically, insert the nozzle into your mouth, align it with teeth or gums, and slightly open your mouth to allow smooth water flow.
5. Press the power button to start. The device automatically shuts off after two minutes. To resume, press the mode button.

Warning

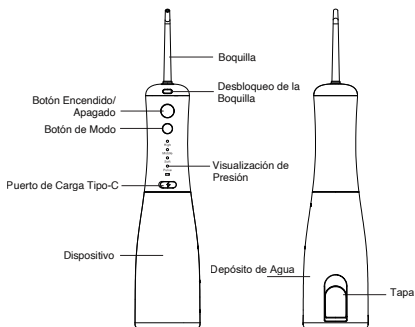1. Only use this product in the modes described in this manual.
2. Do not disassemble the product or replace parts yourself to avoid damage or safety risks.
3. Do not use water hotter than 40°C in the tank to prevent mouth burns.
4. If the product falls into water while charging, unplug the power cord first before removing it from the water.
5. If the plug or power cable is damaged, replace it before reusing the product.

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