
INSTRUCTION MANUAL

Electric toothbrush

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
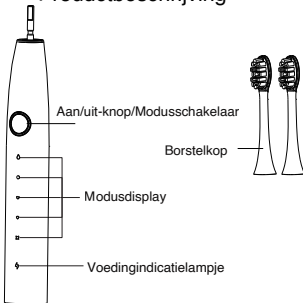
Product Description

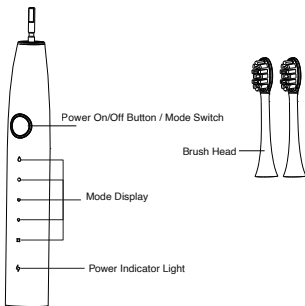Function Description

1. This product features a one-button start system, including on/off and gear switching.
2. 2-minute intelligent timer: after 2 minutes of use in any mode, the product automatically turns off. Restart normally to continue.
3. 30-second zone change reminder: each mode pauses 0.3s every 30 seconds. Press stop or restart to reset the timer.
4. In any mode: press the button once after 5 seconds of operation to shut down.
5. Press the button within 5 seconds to switch modes; press after 5 seconds to turn off.
6. No battery replacement required.
7. IPX7 waterproof rating.

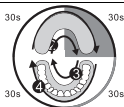
Commutative Diagram

Mode Introduction:

- This product features an intelligent memory function that automatically activates the previous usage mode upon next startup.
- Each mode operates for 2 minutes.
 1. Clean mode - High-speed oscillating bristles easily remove tartar and destroy hidden bacteria on the tooth surface and everywhere (2 minutes).
 2. Sensitive mode - Gentle frequency targeting sensitive teeth and gums (2 minutes).
 3. Massage mode - Very gentle frequency for stimulating and massaging teeth and gums (2 minutes).
 4. White mode - Continuous switching between cleaning and polishing modes helps whiten and maintain the tooth surface, making teeth as white as porcelain (2 minutes).
 5. Polish mode - Ultra-high-frequency power oscillation to polish and brighten teeth (2 minutes).

Preparation for use

1. Align the front part of the brush head with the front part of the brush handle.
2. Attach the brush head to the metal shaft (a small gap remains between both components to prevent the head from colliding with the main unit during operation, enabling the motor to work more powerfully)

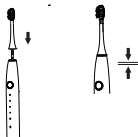
Electric Toothbrush Manual

1. Step one

- Moisten the bristles and apply toothpaste¹

2. Step two

- Press the on/off button to activate the toothbrush.
- On/Off mode switch: In special cases, hold for longer than 5 seconds to turn off. · Apply light pressure during brushing for optimal cleaning.
- Move the brush slowly at a 45° angle between bristles and gums.
- During brushing, you can freely switch between modes via the button.
- The recommended brushing time is 2 minutes.

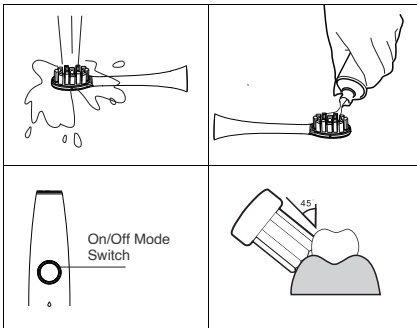Clean the brush head

After brushing, clean the brush head with water (ensure the charging part at the bottom remains covered during cleaning).

Charging


1. When the main battery drops below 3.3V, the charging indicator flashes to indicate charging is in progress.
2. A red light warns of low battery; a green light appears when fully charged.
3. Full charging takes 3-4 hours.
4. If the product has not been charged for a long time, the indicator light may not illuminate immediately during charging. This is normal and usually requires an initial 30-minute charging period before the light functions properly.

Important Reminders

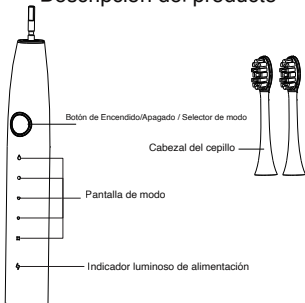1. Do not share the brush head with others.
2. Do not wrap the power cord around the adapter during storage.
3. Avoid placing metal objects (coins, paper clips) on the adapter.
4. Unplug the power cord and store it safely during prolonged non-use to prevent fire hazards.
5. Do not allow children to use this product as a toy.
6. Keep the charging port dry while charging.

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