NEXXT Waterflosser

CO080925

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NEXXT

NL

EN

DE

FR

ES

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FR

ES

Disclaimer: Due to variations in product appearance, if there is a discrepancy between the pixel appearance in this manual and the actual product, please refer to the actual product.

PRODUCT UTILITY

- The strong pulsed water flow can effectively clean bacteria, tooth stains, tartar, and food residues from all parts of the oral cavity.
- The ultra-fine water stream can clean the gingival sulcus and deep between the teeth without barriers, and improve the health of the gums.
- It can effectively remove most of the stubborn plaque on the tooth surface, and the effect is much better than traditional dental floss.
- Quickly clean the oral cavity: Whether at home or eating out, it is the most effective way to quickly restore a clean oral cavity.
- It has a relatively good effect on the peripheral cleaning of orthodontic teeth, dental implants, and dental veneers.
- Help develop healthy daily interdental cleaning habits and keep teeth fresh and clean at all times.

PRECAUTIONS

DANGER

- Please read this manual carefully before use to avoid injury or damage caused by improper operation.
- Do not disassemble this product privately to avoid product damage or electric shock.
- Do not use water above 40 degrees Celsius to fill the water tank to flush teeth to avoid burns.
- If the product falls into the water accidentally during charging, please unplug the power plug before removing the product, and wait until the water stains are completely dry before charging.
- If the power cord or plug is damaged, please replace the power adapter before use.
- Avoid placing the product in a high-temperature place; do not throw the product into a fire, as the battery may explode.

 WARNING

- Do not use the product for any purposes other than oral cleaning.
- Do not use the nozzle to press the gums or teeth too hard to avoid damaging the gums.
- Do not immerse the product in water to avoid product failure.
- Do not use the product in a high-temperature and humid environment.
- Do not continue to use the power cord if it is damaged.


 NOTE

- Do not use this product without the nozzle installed.
- Some functional mouthwashes may be corrosive to a certain extent.
- Please confirm the ingredients of the mouthwash before use, and do not use it indiscriminately.
- If children, the elderly, or those who do not know how to use this product are involved, please make sure that someone supervises and instructs the user.
- Please keep it out of the reach of children and babies, and do not use it as a toy.
- Patients with periodontal disease or those who have had oral surgery in the past 2 months, or who have suffered from oral problems, please consult a doctor first and confirm that it is correct before use.
- Do not use strong gear to continuously wash damaged gums or bleeding teeth to prevent inflammation.
- Please adjust the flushing intensity of the pulse water flow in real time according to your own oral situation.
- Do not soak the whole machine in water; it may cause product failure or complete damage.
- If you do not use this product for a long time, please make sure that the water tank is dry and free of water before storing it to prevent the growth of bacteria.

PRODUCT PARAMETER

Product Weight	: 225g
Water Tank Volume	: 220ml
Charging Voltage	: 5V ---
Charging Current	: 1A
Rated Voltage	: 3.7V ---
Rated Power	: 3.5W

PRODUCT STRUCTURE

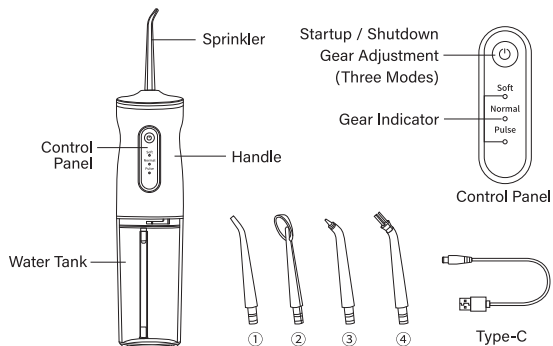
1. The standard nozzle is suitable for daily cleaning.
2. The tongue scraping nozzle is dedicated to cleaning tongue coating.
3. The periodontal pocket nozzle is special for periodontal pocket cleaning.
4. Orthodontic nozzles are suitable for cleaning braces.


OPERATION

- Short-press to turn on. The default mode is strong, and the corresponding indicator lights up.
- Long-press for 2 seconds to switch mode. The mode switches to soft gear, and the indicator light also switches to the corresponding indicator light.
- Switch the mode to the pulse gear, and the indicator light also switches to the corresponding indicator light, and cycles in turn.
- Short-press to shut down. (The instrument will automatically shut down after 3 minutes of operation.)
- When the battery of the instrument is charging, the 3 indicator lights will flash once, and when the battery is full, the 3 indicator lights will be on.


STEPS FOR USAGE

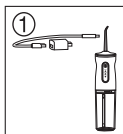
Please fully charge the device before using it for the first time.

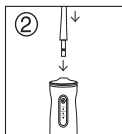
Insert the corresponding nozzle on the handle to ensure firmness and avoid water leakage during use.

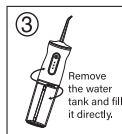Rotate and remove the water tank from the fuselage.

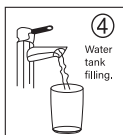Add clean water or mouthwash to the water tank according to your own consumption.

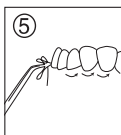Put the nozzle in the mouth and aim at the position to be rinsed.

After use, open the back cover to remove residual moisture.

THE RIGHT WAY TO CLEAN YOUR TEETH


- Choose the gear that suits you. It is recommended to start the trial from the lowest gear and gradually find the gear that suits you best.
- When cleaning the teeth, extend the nozzle into the mouth, aim at the gap between the teeth or the gums, and open the mouth slightly so that the water can flow out smoothly when flushing the teeth.
- After choosing the gear that suits you, align the gums and rinse them in the same direction from the innermost side (such as 1 or 4), until the gums of the last tooth are rinsed, and then turn the nozzle over, align the gums, and continue in the same direction from the inside.
- Rinse in the direction (such as 2 or 3) until it reaches the gum of the last tooth.
- The flushing sequence for the other row of teeth is the same as above.

✓ Demonstrate correctly

✗ Wrong demonstration

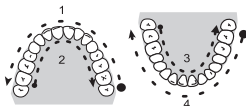WARM REMINDER FOR PRODUCT USE

- Dental experts' suggestion: Use the gentle mode for one week at first use. Then turn to normal mode afterwards.

ATTENTION:

- The direction of water flow is perpendicular to the teeth or gums.
- Move slowly along the teeth and clean them in order.
- The nozzle is aligned perpendicular to the gums.
- Thoroughly clean teeth, trays, crowns, and dentures.
- For the first use, do not spray water directly along the periodontal pocket for a long time, so as not to cause gum discomfort.

TEETH-FLUSHING TIPS

Put the nozzle in the mouth before turning on the switch so that the water will not splash everywhere. In use, the gears are gradually increased from weak to strong to give the gums a stage of adaptation and find a suitable gear that they can bear.

TROUBLE SHOOTING

Common Problem	Possible Causes	Solution
Cannot work or cannot boot.	New purchase or not used for more than three months.	Please charge for more than 4 hours before use.
The water pressure is very low during work.	There is something blocked in the nozzle, or the nozzle is deformed.	Replace the nozzle.
It can only be operated for a short time after charging. The charging indicator does not light up or jump randomly when charging.	The battery has reached the end of its useful life.	Please contact the customer service.
	Insufficient charging time.	Please charge for more than 4 hours before use.
Do not spray water. Cannot work or cannot boot.	The power is not connected.	Please connect on the power correctly.
	Adapter (charging head) does not match.	Replace the adapter or plug directly into the USB socket.
	The product circuit is faulty.	Please contact customer service.
The water pressure is very low during work.	Newly purchased or left unused for a period of time.	Pull out the nozzle to rinse the nozzle hole on the main body.
		Suck hard against the nozzle.

DAILY MAINTENANCE

- If you need to clean the product, please wipe it with water or a damp cloth, and then use a dry cloth or other dry paper towels to wipe off the water left on the surface.
- It is strictly forbidden to soak the product body in water.
- Do not clean this product with hot water exceeding 50 degrees Celsius.
- Do not bend, stretch, or twist the nozzle hose.
- It is recommended to replace the nozzle every 3 months.
- If you do not use this product for a long time, you need to empty the water in the water tank and place it in a cool and dry place.

Correct Disposal of this product

This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

NEXXT

is the innovative brand with high-quality gadgets for in and around the house. Opt for sustainable comfort with our inspiring articles and make your life easier, simpler, or more relaxed. Always at a competitive price. Because we believe in affordable innovation for modern living.

NEXXT, the new ease.

Thank you for choosing one of our NEXXT products!

Before you get started, we strongly recommend reading this entire manual carefully to prevent possible physical injury and damage to the product. We cannot be held responsible for the consequences of incorrect or faulty installation if the instructions in this manual are not followed meticulously. If you have any questions, encounter installation problems, or have doubts about safety, please feel free to contact us at info@conceptonlinebv.nl.

Importer:
Concept Online BV
 De Kreek 28
 NL -3871 MM Hoevelaken
info@conceptonlinebv.nl

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