virtufit

Eye Massager
VF11011

📷 @virtufit_fitness

▶️ VirtuFit

🌐 www.virtufit.com

USER MANUAL

virtufit

VIRTUFIT

Twekkelerweg 263
7553 LZ Hengelo
The Netherlands
info@virtufit.nl

SAFETY INSTRUCTIONS

! WARNING!

Please consult your doctor before starting exercise. This is particularly important for persons with health problems. Read all instructions before using the device. VirtuFit takes no responsibility for any injuries or property damage suffered as a result of using this equipment. To reduce the risk of electric shock, fire or personal injury, this equipment should be used under the following instructions.

- For adult use only.
- This product will block the view. Please do not use it when walking or driving.
- This product has heating function, which must be paid attention to when used by insensitive people.
- Do not use the massage site with skin injury, redness, inflammation, congestion and other discomforts.
- Do not use it for patients with glaucoma, cataract, retinal detachment, eye diseases, high intraocular pressure and physical discomfort.
- After use, the eyes are blurred at first. Close your eyes for 1-3 minutes, and the eyes will gradually become clear.
- Please use the matching USB cable to charge the product.
- Please stop using if you feel unwell during use.

- Do not use it in humid environment such as bathroom.
- Please keep this product away from water, fire and corrosive environment.
- Pregnant women and menstrual women are not allowed to use massagers.
- Before using, please take off the contact lenses, hairpins, hair accessories, etc. to avoid discomfort during massage.
- This product is designed as a health product and is not suitable for diagnosis, treatment and other commercial purposes.
- Please turn off the power switch and take out the battery before discarding the product.
- Please use the product according to the instructions. VirtuFit will not be responsible for any loss or injury caused by improper use.

! WARNING!

Do not use the device without your doctor's approval in the following cases: pregnancy, diabetes with complications such as neuropathy or retinal damage, pacemaker wear, recent surgery, epilepsy or migraines, hernias, recent joint replacements or IUDs, metal pins or plates or other concerns about your physical health. Vulnerable persons and children should be accompanied by an adult when using the device. Research is ongoing on the effects of vibration exercises for specific medical conditions. Vibration exercises in a treatment plan is advisable in some cases. This should be done on the advice of and under the supervision of a doctor, specialist or physiotherapist.


GUARANTEE

The warranty period for this item is 2 years from the day of purchase. The guarantee is only valid when this article is used for activities and normal recreational use and does not cover articles used for rental purposes.

Warranty claims are excluded if the cause of the defect is due to:

- If an item or part is deemed defective by VirtuFit, the warranty covers the repair or replacement of only the defective item or part. VirtuFit is not responsible for any costs, losses or damages incurred as a result of loss or use of this product.
- This warranty does not cover damage caused by misuse, abuse, neglect, normal wear and tear, damage caused by excessive exposure to heat or water, damage caused by improper handling and storage, damage caused by transportation, defective device caused by dropping and damage caused by anything other than defects in materials and workmanship.
- Damage exceeding the cost of the product.
- This warranty is invalid if any unauthorised repair, alteration or modification has taken place to any part of the device.
- The warranty for repaired or replacement parts is valid from the date of original purchase.
- In case of improper use, neglect and/or poor maintenance.
- Failure to maintain the item in accordance with factory instructions.

INSTRUCTIONS

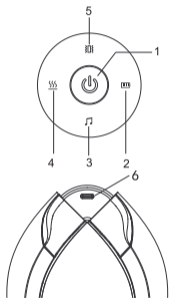1. ON/OFF switch

- Press and hold to turn on and off.
- Short press: Switch Mode.
- Double-click: Turn on and off the music.

2. Battery

- Fast flashing: Power insufficient.
- While charging: Flashing.
- Fully charged: Light stays on.

3. Music

4. Heating

5. Vibration

6. Charging (USB-C Port)

NOTE! Please remove the product protection film before use.

USE

BLUETOOTH

After starting up, the music indicator lights up with a white light. At this time, the Bluetooth of the massager is waiting for connection. The Bluetooth of the mobile phone can be used to search Eye Massager for connection. When the music indicator lights up with a blue light, the Bluetooth connection is successful. When connected to Bluetooth, you can choose your favorite music and adjust the volume.

CHARGING

In the charging mode, whether it is in the startup or shutdown mode, the machine will automatically start the protection program and stop running automatically.

NOTE! Always unplug the charger when the massage head is fully charged.

CLEANING

Before cleaning, turn off the power supply and unplug the power adapter from the socket. Please use a dry soft cloth to clean the surface of the product.

⚠ WARNING! When cleaning, be sure to cut off the power supply. Do not use corrosive detergents such as benzene or diluent to clean products, otherwise it will cause fading, surface aging, structural damage, etc. Use neutral clean sterilized cotton to wipe, and do not use water to clean.

STORAGE

It should not be placed close to the ground, wall or roof. The room should be well ventilated to avoid strong sunlight and corrosive gas erosion.

POWER CHECK

Check the power supply. If the power supply voltage exceeds the specified range, it is not allowed to start up.

TROUBLESHOOTING

Device won't turn on: The battery is low. Please charge it as soon as possible.

No music playing: The music function is off. Double-click the switch button to turn it on.

Flashing charge light: The battery is low. Please charge it as soon as possible.

Continuous "beep" sound: The battery is low. Please charge it as soon as possible.

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