(EN) USER MANUAL: REBBLO - MASSAGE GUN (RB-PC-MSG01-B)

TABLE OF CONTENTS:

- 1) Package contents
- 2) Safety instructions
- 3) Using the device
- 4) Features and specifications

1) PACKAGE CONTENT

- Massage gun
- Massage heads (4 pcs)
- Charging cable
- User manual

2) SAFETY INSTRUCTIONS

Always take the following precautions into consideration before using the "Massage gun" hereinafter referred to as the 'device'.

EXPLANATION OF USED ICONS

	REMARK! Follow the instructions in the user manual!
	WARNING! Warning for potential danger to life and/or serious irreversible injury!
	WARNING! Risk of electric shock Warning for potential risk of death and/or serious irreversible injury!

	WARNING! Risk of suffocation Warning for potential risk of death and/or serious irreversible injury!
	WARNING! Risk of electric shock This device is NOT intended to be used in humid indoor locations. Warning for potential risk of death and/or serious irreversible injury!
	CE MARKING Indicates that the product sold in Europe has been assessed to meet high safety, health, and environmental protection requirements.
	INSTRUCTIONS FOR PROTECTING THE ENVIRONMENT (WEEE, the waste electrical and electronic equipment directive)

IMPORTANT GENERAL INFORMATION

- Read the instructions for use carefully and keep them safe. If you give this device to other people, please also pass on this manual.
- Only use this device in accordance with the instructions in this user manual.
- This device can be used by children under the age of 8 and people with restricted physical, sensory or mental capacity or a lack of experience and expertise, provided they are under supervision or have been given instructions about the safe use of the device and are aware of the potential hazards. Children must not play with the device. Cleaning and maintenance must not be done by children unless they are over the age of 8 and under supervision.

- This device is intended for domestic or similar use, not for professional use.

WARNING! Do not leave packaging material lying around carelessly. This may become dangerous playing material for children. **Risk of suffocation!**

GENERAL SAFETY INFORMATION AND INSTRUCTIONS

- Only use the device for its intended purpose.
- This device is exclusively indented to massage humans.
- Improper use of the device could damage the product.
- The device is not a toy. Keep it out of the reach of children and pets.
- Keep children and pets under strict supervision to prevent them from using the device as a toy.
- Do not allow unsupervised children to clean or maintain this device.
- Never allow children to use domestic appliances as toys or without supervision. Children are unable to judge the hazards associated with incorrectly using electric appliances.

WARNING! Prevent foreign objects, body parts (such as fingers, hair, etc) or garments from interfering with or getting stuck within the device.

- Prior to use, check the device for visible external damage. Never operate the device if it is damaged.
- Do not place the device under any mechanical stress.
- Please handle the device carefully. Jolts, impacts or falls, even from a low height, can damage the device.

- Do not use or store the device in an environment subject to strong static electricity or magnetic fields.
- Do not use the device in the vicinity of flammable materials.
- Do not use the device in the immediate vicinity of heat sources such as naked flames (hobs, for example).
- Do not expose the device to direct sunlight, heat sources, excessive humidity or corrosive environments.

WARNING! This device is **NOT** intended to be used in humid indoor locations.

- The warranty will expire immediately if the device is damaged as a result of not observing this user manual. The manufacturer accepts no liability for damage caused by non-compliance with the user manual, negligent use or use which does not comply with the conditions of this user manual.
- Do not use abrasive cleaning solvents to clean the device.
- Do not use foreign objects to clean the device.
- Do not use the device if there are visible signs of damage.
- Do not use the device if it has been dropped or does not operate normally.
- Do not attempt to disassemble or repair the device or modify it in any manner.
- The device does not contain any parts which can be serviced by the user.
- Only allow repairs to be carried out by a qualified service centre. Qualified service centre: after-sales department of the manufacturer or store, which is recognized and authorized to carry out such repairs, so that potential hazards can be prevented. If

problems are encountered with the device, please return it to this department.

BATTERIES, ELECTRICITY & HEAT

- The device's rechargeable battery should be charged for ± 5 hours before using the device for the first time.
- Make sure your hands are dry before touching the device or charging cable.
- When plugging in the charging cable, make sure people cannot accidentally pull it out or trip over it.
- Never leave the device to charge unattended.
- Regularly check the charging cable for damage.
- Do not charge the device if the charging cable shows signs of damage.
- Make sure the charging cable cannot be trodden on or pinched.
- Do not allow the charging cable to hang over sharp edges and keep it away from hot objects and naked flames.
- Remove the charging cable from the device and the power adapter from the electrical outlet during lightning and storms.
- Make sure that the device and charging cable do not come into contact with heat sources, such as hotplates or naked flames.
- Once charging has been completed, always disconnect the charging cable from the device.
- Do not attempt to charge the device when the battery is damaged.
- To reduce or prevent the risk of personal injury, electrical shock, fire or damage to the device, never insert metal objects in any opening of the USB type-C charging port or

any other openings or ports.

- The rechargeable batteries are built into the device and cannot be replaced.

WARNING! Do not open the casing as this may result in electric shock!

WARNING! To avoid the risk of electric shocks, never immerse the device in water or any other liquid. And do not use it again if this happens.

- Do not crush, puncture or dismantle the device; this might damage the battery.

WARNING! Damaging the casing and/or piercing the device, can lead to explosion or fire!

WARNING! To reduce the risk of fire or electric shock, protect the device from dripping or splashing water and keep away from objects filled with liquids, such as vases. Such objects should not be placed on or near the device.

- Do not use the device in temperatures that fall below -10°C or exceed 40°C .
- Battery cells can leak under extreme conditions. Wear protective gloves when handling leaked batteries and dispose of the gloves after use. If the battery leaks, do not allow any battery fluid to come into contact with your skin or eyes. If contact has been made, wash the affected area with copious amounts of water and seek medical advice.

WARNING! Do not throw the device into a fire, as it may explode!

- Do not dispose of the device in your household waste. Take it to your local recycling station. Contact your local government offices for details.

WARNING! Never cover the power adapter, charging cable or device with cloths, curtains or any other material in

order to prevent overheating and potential fire hazards.

IMPORTANT INFORMATION, FOLLOW THESE INSTRUCTIONS

- Make sure only one person uses the device at a time.
- Do not use the device for more than 10 minutes continuously.
- Do not use this device if you are affected by one or more of the following conditions and illnesses: circulation disorders, varicose veins, inflammation of the veins or thrombosis, open wounds, skin sensitivity (such as diabetes, eczema, psoriasis), skin irritation or rashes, bruises, broken skin, reduced pain sensitivity, or if you are taking painkiller or tranquillizer tablets.
- Do not use the device if you are tired, ill or under the influence of alcohol, drugs or medication.
- If you experience any pain, discomfort or swelling using the device, or find the massage unpleasant, stop the treatment and consult your doctor.

WARNING! Consult your doctor in all cases before using the device under the following circumstances: a heart disease, pacemaker, cancer, intra-articular fixation, artificial joints and/or (electronic) implants, aneurism, blood thinners.

- Never massage the heart area with the device if you have a pacemaker or any other electrical medical device.
- Do not use the device on any particularly sensitive part of the body, such as the head,

face, larynx, spine or any hard or bony area of the body.

- Do not use the device if you suffer from a herniated disc or from abnormal changes to the vertebrae.
- Only use the provided accessories with the device.
- Only use on a dry clean surface of the body by lightly pressing and moving across the body for approximately 60 seconds per region.
- Bruising can occur regardless of speed levels or pressure applied. Check treated areas frequently and immediately stop at first sign of pain or discomfort.
- Do not use the device during pregnancy.
- Do not use the device when you are in bed or sleeping.
- Do not use the device in a driving vehicle.
- Under no circumstances use the device in combination with (essential/massage) oils, lubricants or other similar substances.
- Not for internal use, for external use only.


3) USING THE DEVICE

Open the packaging carefully before taking out the device. Make sure the device is complete and undamaged. If any components are missing or damaged, contact the vendor and do not use the product. Retain the packaging or recycle it in accordance with local regulations.

Before use, it is advised to clean the device thoroughly with a damp cloth and dry it thoroughly with a dry cloth.

PARTS & FEATURES OF THE DEVICE

- A. Massage heads (4 pcs)
 - A1. Round massage head
 - A2. Flat massage head
 - A3. Bullet massage head
 - A4. Fork massage head
- B. Massage head holder
- C. LED indicator
- D. Power button
- E. Charging port
- F. Handle
- G. Charging cable

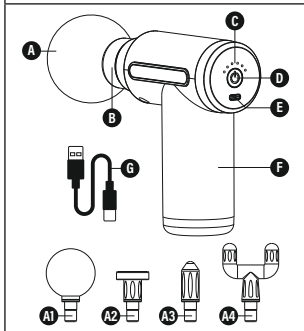- indicator lights (C) will light up and the device will automatically stop charging.
- Unplug the charging cable from the device, which is now ready for use.
- To check the battery level, press the power button (D) once when the device is powered off and the LED indicator (C) will light up.
- The charging LED indicator (C) has 6 lights. Each light represents 1/6 of the battery.
- If the battery level is low, only the first light of the charging LED indicator (C) will light up, indicating the battery needs to be recharged.

INSTALLING THE MASSAGE HEAD

To use the device, it must be fitted with a massage head. Different muscle groups require different care. Select the desired massage head.

- Massage heads (A1) & (A2) are suitable for large muscle groups and body parts.
- Massage head (A3) is suitable for precise pain points.
- Massage head (A4) is suitable for specific muscle groups such as next to the neck, spine, waist and Achilles tendon.

Push the selected massage head (A1), (A2), (A3) or (A4) on the massage head holder (B). To change or remove the massage head, make sure the device is powered OFF and gently pull the massage head off the massage head holder (B).

CHARGING THE DEVICE

- Make sure the device is turned off.
- Carefully plug the connector of the charging cable (G) into the charging port (E).
- Plug the USB connector of the charging cable (G) into the USB port of your USB charger.
- The LED indicator (C) flashes when the device is charging.
- Once charging is complete, all the LED

Note! Never use the device without a massage head attached.

USING THE DEVICE

Before use make sure you follow and respect the instructions as described in "Important information, follow these instructions". Power on the device by pressing the power button (D) for ± 3 seconds. Once the device is

on, press the power button **(D)** to cycle through the 6 different speed levels, from low to high. Hold the device firmly by the handle **(F)** and place the massage head **(A)** on the part of your body you want to massage. Move the device around accordingly. **Use the device a maximum of 60 seconds per region and observe a total operation time of 10 minutes.** Power off the device by pressing the power button **(D)** for ± 3 seconds.

Note! Do not use the device whilst charging.

CLEANING THE DEVICE

Proper and regular cleaning will ensure safe operation and extend the lifespan of the device.

- Make sure the device is powered off and the charging cable is unplugged.
- Use a damp cloth to remove dirt from the device and dry with a dry cloth.
- Do not clean the device with strong chemicals - or alkali-based, abrasive or disinfecting agents - because they could damage the surface of the product.
- Do not immerse the device into water or any other liquids.
- Store in a dry place.

4) FEATURES AND SPECIFICATIONS:

- Input: 5.0 Vdc, 1.0 A (max.)
- Battery: 1.800 mAh, Li-ion, 7.4 Vdc (13.32 Wh)
- Power: 24 W
- RPM: 1.800 – 3.000 RPM
- Weight: 485 g
- Dimensions: 14 x 5.5 x 14 cm (H x W x D) without massage head

All contents and specifications mentioned in this manual are subject to change without prior notice.

INSTRUCTIONS FOR PROTECTING THE ENVIRONMENT

(WEEE, the waste electrical and electronic equipment directive, applicable in EU member states)

Your product has been manufactured using high quality materials and components which can be recycled and used again. At the end of its service life, this product cannot be disposed of in normal domestic waste, but must be taken to a special collection point for recycling electrical and electronic equipment. This is indicated in the user manual and in the packaging by the crossed-out wheeled dust bin symbol. Used raw materials are suitable for recycling. You can ask your local competent council about the location of your nearest collection point. By recycling used appliances or raw materials, you are playing an important role in protecting our environment.

Correct disposal ensures that waste electrical and electronic equipment (which may contain hazardous substances) is recycled and reused appropriately. It helps avoid potential damage to the environment and human health, wildlife and to preserve natural resources.

Batteries:

Do not dispose batteries in your household waste. Dispose of batteries by taking them to your local recycling station. Contact your local competent government offices for details.

Following information is only valid for EU member states.

Pb, Hg, Cd If the battery contains higher values of lead (Pb), mercury (Hg) and/or cadmium (Cd) than defined in the Battery Directive (2006/66/EC), the chemical symbols Pb, Hg and/or Cd are shown below the crossed-out wheeled dust bin symbol.

EU DECLARATION OF CONFORMITY

We,
Mascot Online B.V.
Veluwezoom 32
1327 AH Almere
The Netherlands

Declare under our sole responsibility that the product:

- Product Brand: Rebblo
- Product Model No: RB-PC-MSG01-B
- Product Description: Massage gun
- Product Rating: Input: 5.0 V, 1.0 A max.
- Battery: 1.800 mAh Li-Ion (13.32 Wh)

Following the provisions of the Directives:

- Electromagnetic Compatibility (EMC): 2014/30/EU
- Restriction of the use of certain Hazardous Substances in Electrical and Electronic Equipment (RoHS) Directive: 2011/65/EU and amendment directives (EU) 2015/863

And standards to which Conformity is declared:

EMC:

- EN IEC 55014-1:2021
- EN IEC 55014-2:2021
- EN IEC 61000-3-2:2019+A1:2021
- EN 61000-3-3:2013+A2:2021

Place: Almere
Date of issue: 01-08-2022
Mascot Online B.V.

Authorized Signature name: B. van Dijk
Title: Managing Director

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