
ITEM: WH-WMC

OWNER'S MANUAL

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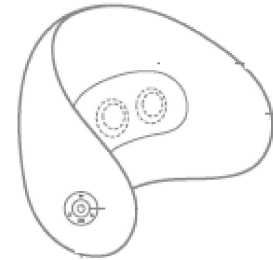
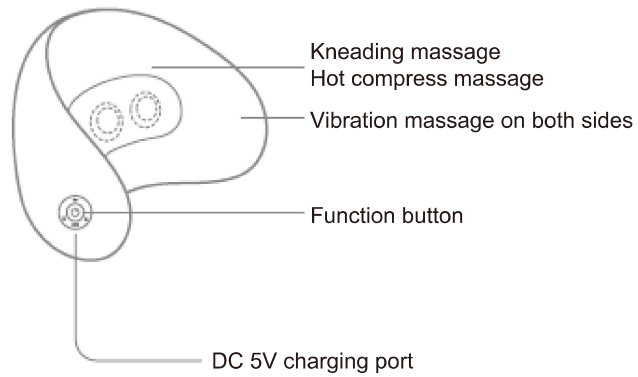
Note: Please use a charger with a rated voltage of 5V for charging. It is recommended to use a 5V/2A charger.

On/Off
Pressing the On/Off button will turn on the kneading and vibration massage functions at the same time. It automatically stops working after 15 minutes.

Hot compress massage
Press to turn on the hot compress. Press again to turn off the hot compress.

Kneading massage
Switch two modes. Mode 1 (default) is strong, mode 2 is gentle.

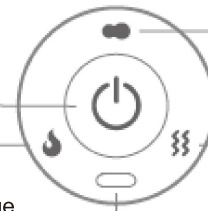DC 5V charging port


Vibration massage
Switch two modes. Mode 1 (default) is rhythm, mode 2 is continuous.

CAUTION


1. In order to prevent skin abrasion, do not use the product in direct contact with the skin.
2. If you feel uncomfortable using this product, please stop using it immediately. If you feel seriously uncomfortable, please seek medical advice in time.
3. It is a home health product, which is not suitable for diagnosis, treatment and other medical uses.
4. Poor contact of the charger may cause the product to stop working.
5. Improper use of the product may damage your body and cause accidents.

Attention

1. You must use a charger with a rated voltage of 5V.
2. Do not pull the charging cable and charger connected to the product too hard.
3. Do not use it in a bathroom or damp places.
4. To avoid product damage, excessive noise during working and prevent a decline in durability of the product, please do not use the product on solid objects.
5. The product has a heating surface. People who are insensitive to heat must pay attention during using.
6. If you find any abnormal about the product, please do not disassemble it by yourself, and contact the customer service.

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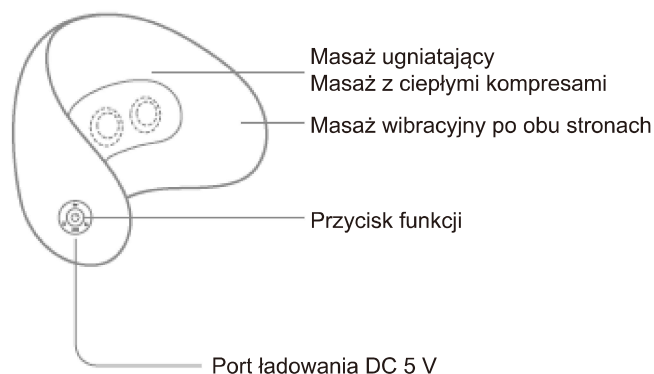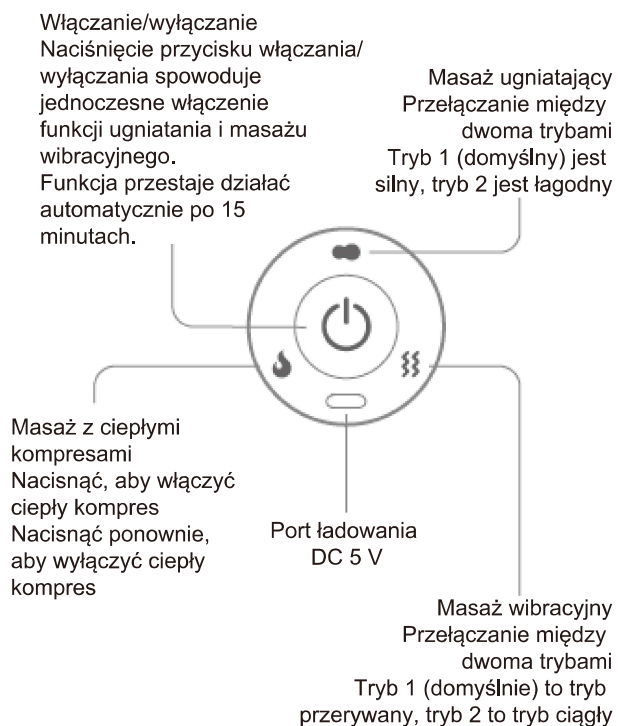DECLARATION OF CONFORMITY

Confinity NV, Dorp 16, 9830, Sint-Martens-Latem, Belgium
declares the following device in sole responsibility:

Brand name: Wellhealth

Product type: Wellhealth waist massager comfort

Item number: WH-WMC

Complies with the following harmonization rules:

EN IEC 55014-1:2021, EN IEC 55014-2:2021

EN IEC 61000-3-2:2019/A1:2021,

EN 61000-3-3:2013/A1:2019

Signed for and on behalf of:

Sint-Martens-Latem, Belgium - MAY 2023

A. Pappijn - Product Manager

A handwritten signature in black ink, appearing to read 'A. Pappijn', is written over a horizontal line.