WellHealth

INSTRUCTION MANUAL

PREMIUM ORAL CLEANER

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The device can be flushed with water	can't be soaked in water for a long time


1. Product Description

This device is a real ultrasonic tooth cleaner that uses a high-frequency ultrasonic transducer to generate high-frequency vibration. After turning it on, you won't feel any vibration when holding it in your hand. Only when it touches hard objects (such as teeth). You will feel the vibration (it won't work when it touches your gums). The operation is simple, convenient, fast, and efficient and can effectively remove dental calculus, tobacco stains, plaque, and tea stains.

2. Attention

1. After the device is turned on, there is no sound except the indicator light will light up. When the cleaning head touches the teeth, it will make a beeping sound, which is normal operation.
2. When the cleaning head touches the boundary of the gums, it is normal that gum bleeding may occur even in the case of healthy gums. But if bleeding continues, it may be due to gum defects, such as periodontal disease. Please stop using it and consult your dentist.
3. Do not leave the cleaning head on the surface of the teeth for a long time. Please use the cleaning head to scrape off at 45° to remove stubborn calculus and tooth stains.

3. Packing List

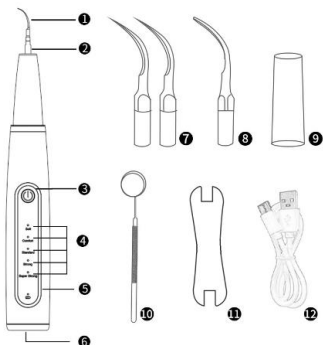① Cleaning Head	② LED Light
③ Power Button	④ Mode indicator
⑤ Power Indicator	⑥ USB Charging Port
⑦ Pointed Cleaning Head	⑧ Flat Cleaning Head
⑨ Dust cover	⑩ Oral Mirror
⑪ Wrench	⑫ Charging cable

4. Specification

Charging Voltage	5V
Charging Current	≤1A
Charging Time	2-3 Hours
Frequency	32KHZ-40KHZ
Consumption	≤2.8W
Battery Capacity	1500mAh
Weight	140g

5. Charge Before USE

Please fully charge the device before use (fully charged in 2-3 hours)

5.1. The power indicator light is red when charging.

5.2. The power indicator light is green after being fully charged.

5.3. During the charging, the device cannot work.

5.4. Fully charge in 2-3 hours. It is recommended to use a 5V/1A adapter.

6. Operation Guide

6.1. Power on the Device:

Press the power button to turn on the device. The indicator light is green, the green light in mode 1 is on. to enter the ultrasonic soft mode. In use, press the Power Button again to enter mode 2 (normal mode), and the green light in mode 2 is always on. By analogy, there are 5 modes in total:

(1)Soft;(2)Comfort;(3)Standard;(4)Strong;(5)Super Strong

6.2. Turn On LED Light:

Turn on the LED Light: After starting the device, press the LED Control Button and the LED light will light up. If you press the light switch again, the LED light will turn off.

6.3. Power off:

After use, press and hold the power button for 2-3 seconds, the indicator light goes out and power off. The device has the overheat protection function enabled, and it will automatically power off after 10 minutes of work. If you want to continue using it, just start it.

6.4. During use, if the power indicator is red, it means that the power is too low. You need to charge the device to continue using it. If the red light flashes, that means the battery is run out. Also needs to be recharged.

7. Introduction

7.1. Pointed Cleaning Head: used to remove calculus and tooth smoke stains on gingiva.


7.2. Flat cleaning Head: used to remove large stones and smoke stains on the neck and lingual side of the supragingival teeth.

⚠ To prevent cross-infection, do not share the cleaning head with many people.

8. How to Use

The cleaning head gently touches the teeth vertically and moves to the top and bottom as well as left and right. Please adjust the use time of each part according to the condition of the teeth.

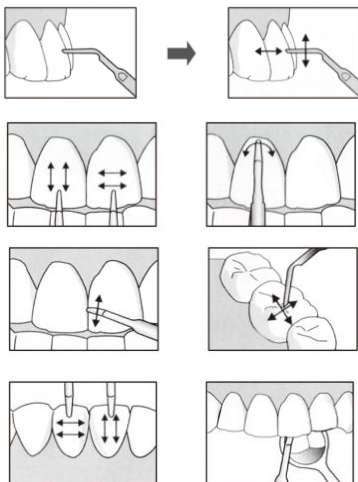
9. How to Install


The orientation of the cleaning head may be different, but as long as it is firmly mounted on it, it can be used in any direction.

Insert the wrench
in position A

confirm position A

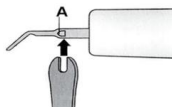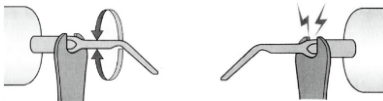Rotate the hend counter Clockwise with a wrench to separate the working tip and complete the disassembly

Avoid over-tightening, which can damage the interface.

10. How to Clean

10.1. Before the cleaning head is used, it can be disinfected by steaming, alcohol soaking, high temperature and high pressure, etc.

10.2. When not in use for a long time. Please keep it in a dry place.

10.3. Please do not use diluent, benzene, or solvents, such as kerosene or scrubbing powder.

The device can be flushed with water	can't be soaked in water for a long time

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