
WellHealth

INSTRUCTION MANUAL

PREMIUM TOOTH POLISHER

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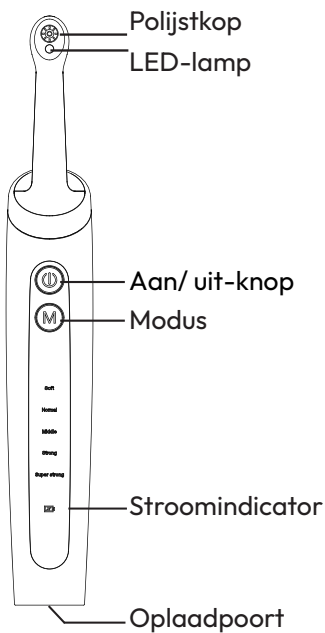2. Extra

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Puntige
borstelkop


Platte borstelkop

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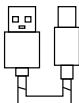Oplaadkabel

Sleutel

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Stap 1

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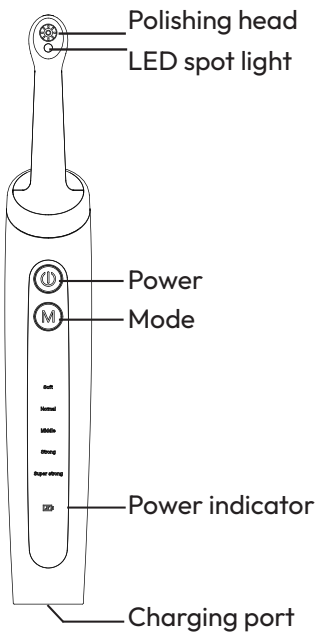
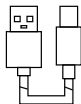
3. How to use

Step 1

Tighten the cleaning head clockwise by hand.

NOTE: Choose the proper teeth cleaning head.

(1) Silicone Rub Head: Apply directly to the tooth surface to remove stubborn stains, plaque, and bacteria.

(2) Flat-top Brush Head: Apply directly to the tooth surface to remove sanding residues.

(3) Polishing Cup Head: For tooth surface polishing with toothpaste/polishing paste (not included).

(4) Umbrella Polishing Cup Head: For interdental, behind tooth polishing with toothpaste/polishing paste (not included).

Step 2

Short press the Power Button to turn it on, during the process, you can press the Mode Button to adjust the speed according to your needs.

Step 3

Silicone Rub Head: mim the tip of the silicone rub head on stubborn stains

NOTE: (1) Keep moving; lift it every 2-3 seconds; the polishing effect is better; the silicone head will not damage your teeth, but it will heat up for a long time.
(2) Some powder may float out during the process, which is a normal phenomenon.

Step 4

Replacing the head

Use a wrench to hold the metal part of the head of the machine, and then unscrew it counterclockwise by hand.

Step 5

Apply directly to the surface of the teeth to remove the residue after grinding; just like brushing your teeth, you can use a small amount of toothpaste in conjunction with.

Step 6

Polishing Cup Head

Squeeze the toothpaste or polishing paste into the cup, turn on the machine, align the cup with the tooth surface, start on the low gear.

(1) If your toothpaste has a thickness, please use the paste from coarse to fine, use up all the coarse ones, then use the next model and repeat the polishing steps for the best effect.

(2) When squeezing the toothpaste, just fill it up; otherwise, it will splash out.

(3) Keep moving; lift it every 1-2 seconds; the polishing effect is better; the silicone head will not damage your teeth, but it will heat up for too long.

(4) Continue to replenish toothpaste in the cup; it is recommended to replenish it every 5 seconds to ensure the whitening effect.

(5) A lot of saliva mixture will be produced in the process; please clean it up in time to prevent splashing.

(6) Polish each tooth in turn, and the effect is better.

The recommended time for each tooth is 10s.

Step 7

Umbrella Polishing Cup Head

In the same way as the Polishing Cup Head, the cleaning position is between the teeth and behind the teeth, and the umbrella shape is more suitable.

Step 8

Short press the switch button to turn it off

NOTE: The device will shut down after 3 minutes of continuous operation. Just restart it for further cleaning.

Step 9

Rinse your mouth with water and spit out

NOTE: We don't recommend smoking or eating dark food like coffee and tea within 7 days after cleaning.


Step 10

Wipe and clean the handle of the device with clean, wet cloths. Then put the cover on and keep it out of the reach of children.

4. Safety

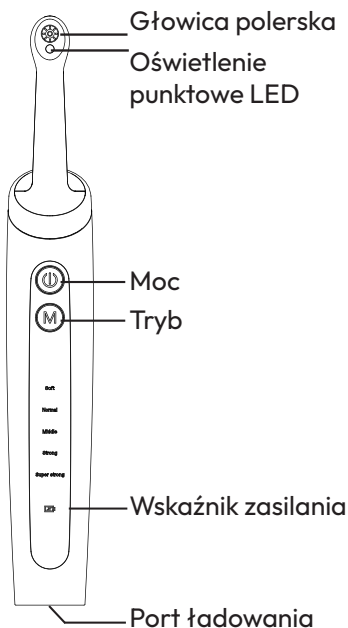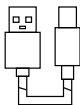
1. This device is not intended for use by children (less than 15 years old).
2. Keep the power cord away from the heated surface. The temperature of the product may rise slightly when it is working or charging.
3. Do not drop the device, or insert any objects into any opening or hose.
4. Do not operate the device with a damaged cord or not work improperly.
5. Only used on the tooth surface to prevent touching the gums and causing unnecessary bleeding.
6. Do not use pressure. This product is engineered to clean/polish appliances as well as teeth; use caution and do not overpolish. A light touch is all that's needed.
7. It is normal that a loud sound will be produced during the use of this product.
8. Never use the device while charging.
9. We recommended starting from the soft mode when first using it.

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